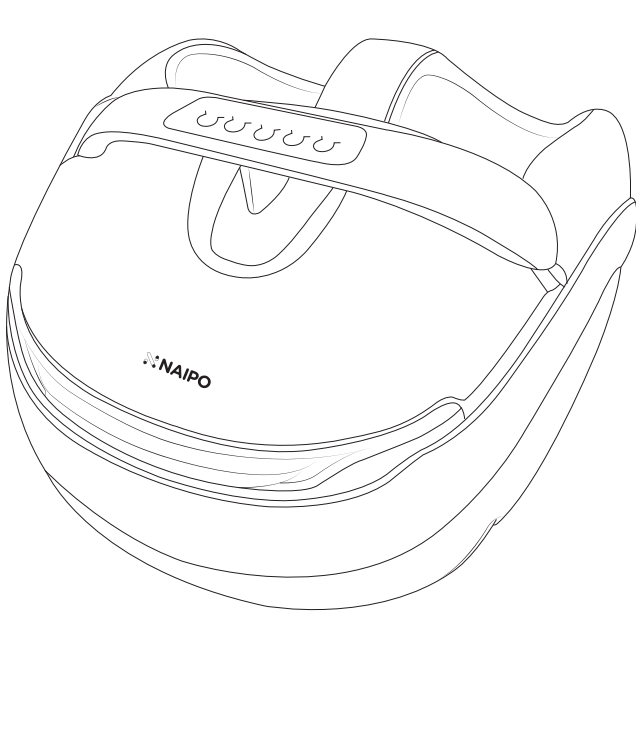
## User Manual

### Foot Massager With Heat and Rolling Massage

Model: MGF-839A

Shenzhen Pinlong E-Commerce Co., Ltd.  
Web: www.naipocare.com  
E-mail: support@naipocare.com  
Add: Huanan City Electronic Trading Building P07  
Room 120, Longgang District, Shenzhen, China.

#### Introduction

Thank you for purchasing the Naipo Foot Massager With Heat and Rolling Massage. Please read this user manual carefully and keep it in a safe place for future reference. If you need any assistance, please contact our support team with your product model number and Amazon order number.

#### Included

- 1 x Foot massager
- 1 x User manual
- 1 x Adapter

#### Product Specifications

|                   |   |
|-------------------|---|
| Name              | Naipo Foot Massager With Heat and Rolling Massage |
| Model             | MGF-839A  |
| Voltage           | 100V-240V 50Hz/60Hz                               |
| Power consumption | 65W   |
| Timer             | 15 minutes  |

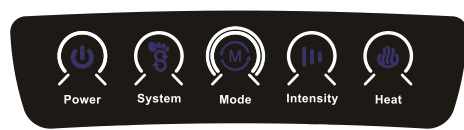
#### Important Safeguards

- Please avoid the direct contact of your skin with the product.
- The product is equipped with heating device, so be aware of this if users are people insensitive to heat.
- Please do not use the product in an environment where there is a high percentage of gasoline, diluent, alcohol, or other corrosive detergents like bleach to avoid color changing, deformation, cracks, or malfunction.
- Please do not sit on the product or place heavy objects on it.
- Please do not expose the product to direct sunlight.
- People listed below should not use the product:
  - Pregnant women;
  - People who suffer from heart diseases or other health problems.
- It is recommended not to use the product for more than 15 minutes at a time. Overly long use of a time might cause muscle pain or nerve pain.
- Children are not allowed to use the product without supervision from adults.
- Before or when using the product, please do not apply medical/cosmetic cream or emulsion on the product, or use any foot powder.
- Please do not use the product if there is damage on foot sleeves or the power cable.
- Please do not use the product if there enters inside it by accident.
- Please stop using the product if abnormal or severe pain occurs.

#### Product Diagram

- Zipper
- Control panel
- Lower shell
- Detachable foot sleeves (detaching: unzip to detach foot sleeves / installing: attach the hook and loop fastener and then zip up)
- Upper shell

#### Control Panel

**Power:** When the product is under the standby mode, the red light flickers; Press the button to turn on the product with functions like automatic massage program, heating function, and tapping massage being activated (green light); Press the button again to stop all functions.

**System:** When the product is working, press this button to shift among 3 automatic massage programs (with green light, blue light, and red light interchanging for each program).

**Mode:** Under automatic massage program, press this button to control tapping massage, rolling massage, and air compression.  
When the product is turned on, tapping massage, air compression, and rolling massage will be activated (green light); Press the button to stop tapping massage while air compression and rolling massage remain (blue light); Press this button again to stop rolling massage while air compression and tapping massage remain (red light); Press the button for the third time stop rolling massage and tapping massage while air compression remains (light: off).

**Intensity:** Press the button to adjust among 3 intensity levels (low: green light / medium: blue light; high: red light); When the product is turned on, low intensity level would be activated by default.  
**Note:** Long press the "Mode" button to activate tapping massage alone; Press the "Intensity" button to adjust tapping massage speed among 3 levels; Press the "Mode" button to choose among 3 tapping massage modes.  
If tapping massage is not activated alone, press the "Intensity" button to adjust air compression among 3 intensity levels.

#### Product Maintenance

- Please disconnect the product from the power supply before cleaning.
- Please clean only with soft and mildly damp sponge.
- Please do not clean with any corrosive detergent (like benzene and diluent).

#### Questions & Answers

| Problem   | Solution   |
|---|--|
|   | Please check whether the power cable is in good condition without any damage.  |
|   | Please check whether the plug is connected properly.   |
| The product fails to work after the power supply. | After working for long, the temperature of the product and the motor might get overly high, and the product will automatically shift to protection state; Please disconnect the product from the power supply and do not use for 30-50 minutes until it cool down. |
| The product stops working during the use.         | If the timer has been activated and time is up, please restart the product.<br>Disconnect the product from the power supply, reconnect it to the power supply, and then restart the product.   |

#### Warranty & Customer Support

For questions, assistance, support, or warranty claims, please contact us via the email or the address below.  
Web: www.naipocare.com  
E-mail: support@naipocare.com  
Address: Huanan City Electronic Trading Building P07 Room 120, Longgang District, Shenzhen, China.

\*Note: Products purchased directly from Naipo benefit from our after-sales service. If the purchases are made from other selling firms, we recommend you contact them for warranty issues.

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Web: www.naipocare.com  
E-mail: support@naipocare.com  
Address: Huanan City Electronic Trading Building P07 Room 120, Longgang District, Shenzhen, China.

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