




User Manual

Foot Spa Massager with Full Roller,
Heat and Bubbles

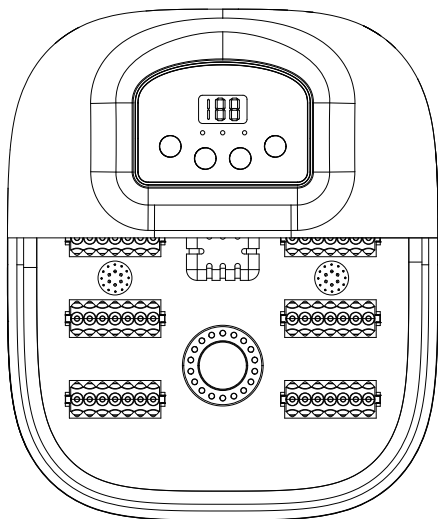

Model: NPFT-SPA01

Contents

English	01~04
---------------	-------

Introduction

Thank you for purchasing the Naipo Foot Spa Massager with Heat, Bubbles and Vibration. Please read this user manual carefully and keep it in a safe place for future reference.


Included:

- 1 x Foot Spa Massager
- 1 x User Manual

Product Specifications

Model	NPFT-SPA01
Name	Foot Spa Massager with Full Roller, Heat and Bubbles
Rated Voltage	220~240V (EU) / 120V (US)
Rated Frequency	50Hz (EU) / 60Hz (US)
Power Rating	500W
Net Weight	1.677 Kg
Product Dimensions	380*330*200mm
Temperature Range	95°F-118°F (35°C- 48°C)

Product Diagram

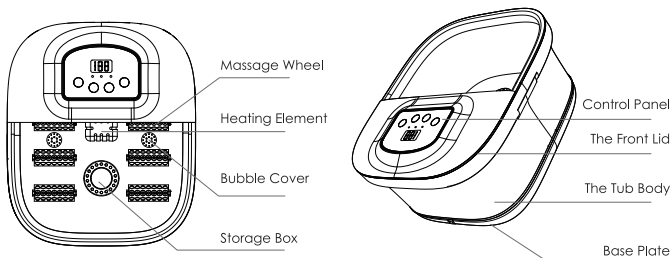

Control Panel



User Tips:

1. The range of water temperature is recommended to be 95~113°F/35-45°C (overly high temperature might hurt the skin).
2. Never turn on the heating function if there is no water in the product, other wise the motor would be burnt, leaving the product unable to work any more.
3. The heating element is located in the center of the product, so it's strongly recommended to turn on the bubbling function while the heating function is on for accelerated circulation of water heating.

Q & A Product Operation Guide

- Lay the product flat and fill it with moderate amount of water (the temperature of water is recommended to be 95°F (35°C) to save your time for waiting). Connect the product to 120V/60Hz (240V/50Hz) power supply and "OF" will appear on the screen.
- Press the power button and the temperature of water is shown on the screen.
- Press the heating button and set your favorable heating temperature. The heating temperature will rise by one degree with each press. When the set temperature is 118°F (48°C), it will reverse back to 95°F (35°C), which is a cyclical program. The default heating temperature is 107°F (42°C) with a range from 95°F-118°F(35-48°C). To heat the water evenly and improve comfort level, it's recommended to turn on the bubbling function at the same time.

- Press the vibrating button to activate vibrating massage and press the button again to stop vibrating massage.
- Press the bubble button to activate bubble massage and press the button again to stop bubble massage.
- You can put foot spa materials such as Ay Tsao, rose petals and salt bags in the storage box. "L/R" is marked on the storage case, so please pay attention to the direction when placing.

Q & A

Problem	Cause	Remedy
The heating function doesn't work.	The product is not connected to the power supply.	Please check and connect the product to the power supply .
The product doesn't work.	The product is not initiated properly though connected to the power supply.	Please use the product strictly following the user manual.
The product stops working during the use.	The product works beyond default time.	Please restart the product.

Important Safeguard

01. Please do not use the product if it is damaged.
02. Read the user manual carefully before the product is put into use.
03. Please do not use the product without water inside the tub and keep the water volume between minimum line and maximum line. Please disconnect the product from mains supply after use.
04. Please do not stand on the product or put heavy things inside it.
05. Please do not pull the power cable and pull or lift the product by the cable, otherwise it might get damaged.
06. Please make sure that the product is always positioned horizontally.
07. Please stop using the product when uncomfortable feelings occur, as feeling dizzy and headache.

08. Please make sure that the cable line and plug are dry and keep them away from fire when the product is in use.
09. Please do not fill the tub with liquid of acidity, basicity or corrosivity and liquid containing particles or floccules.
10. Please do not fill the tub with water above 113°F (45°C), otherwise your feet might get burnt or the product might be out of shape because of high temperature.
11. This product is not recommended for people with inflammation, trauma, ringworm and burns on feet.
12. People with heart disease, people with serious embolism and cerebral hemorrhage, pregnant women and people with terminal diabetes are not recommended to use the product without consulting a doctor.
13. Children, people who are disabled (both physically and mentally), sluggishers and those who are insensitive to heat or cold are recommended to use the product only with assistance from others.
14. Please do not use this product before or after meals within 30 minutes.
15. Do not use without water
16. Do not use while standing

Massage Areas

The foot spa massager is for foot spa and foot massage.

For questions, assistance, support, or warranty claims, please contact us via the email or the address on the cover.

***Note:** If the purchases are made from other selling teams, we recommend you contact them for warranty issues.

Shenzhen Bestvap Technology Co., Ltd.

Web: www.naipocare.com

E-mail: wholesale@naipocare.com

Add: 503 ManTao Building 720 Minzhi Avenue Minzhi Street Longhua District, ShenZhen, GuangDong 518000