




Smart handheld massager MGPC-806P

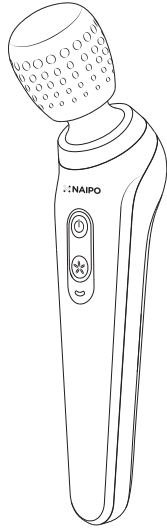

Instruction Manual

Thank you for purchasing our handheld massager. Before use, please read this manual carefully and pay special attention to safety matters to ensure correct use.

Safety Precautions

Mark Definition

	Indicate that you may cause personal injury or property damage due to misuse.
	Indicate that you may damage the machine or may be injured by misuse.
	Prohibit use
	Prohibit disassembly

1. For the following people, please consult your doctors before use.

	<ul style="list-style-type: none"> ● People implanted cardiac pacemaker ● Pregnant or woman shortly after giving birth ● People suffering from malignant tumor ● Do not use this product while charging ● The device should not be disassembled or modified by itself, otherwise it may cause damage or injury ● After charging, disconnect the USB cable, otherwise, long-term power-on will shorten the battery life of the device 	<ul style="list-style-type: none"> ● People suffering from heart disease or strained back disease ● People suffering from fractures, osteoporosis or arthropathy ● People suffering from skin trauma or eczema ● Do not let children use or play directly to avoid injury ● If any abnormality or failure is found during use, please stop using it immediately. ● This product is not waterproof, please use it away from water
--	--	--

2. Things must know when using


	1. This product is not used for medical purposes, but only provides personal massage and health benefits.
	2. Do not allow children or people without autonomy to use this product. Please keep this product out of reach of children.
	3. When the product is plugged in, there must be an adult nearby, especially if there is a child.
	4. Do not use this product when the following situations happen: The data cable damaged, abnormal use, dropped off the ground or into water, damaged.
	5. Do not stand or sit directly on the product, because it may cause fall accidents.
	6. Excessive use will cause the product to overheat and shorten its life, so please cool it down for a while before reusing.

	<ol style="list-style-type: none"> 1. Do not use it in humid environments such as bathroom. 2. Do not immerse it in water or other liquids to avoid short circuit or electric shock. 3. Do not place the product near stoves or direct sunlight. 4. Do not use it during lightning strikes. 5. Do not insert or drop any object into any opening of it.
	<ol style="list-style-type: none"> 1. Do not privately disassemble the internal components of this product to avoid any possible machine failure. 2. The product cannot be disassembled, repaired or modified except our designated repair technicians.

3. Charging safety

	<ol style="list-style-type: none"> 1. Do not connect metal conductors with the battery poles, or touch the terminals of the battery, so as to avoid short-circuiting the battery, as well as burns or other physical injuries caused by the battery overheating. 2. Do not expose the battery to high temperatures around heating devices, such as sunlight, heaters, microwaves, ovens or water heaters, etc. Overheating of the battery may cause explosions. 3. Do not disassemble or modify the battery, insert foreign objects, or immerse in water or other liquids, so as not to cause battery leakage, overheating, fire or explosion. 4. If the battery is discolored, deformed, or abnormally overheats during use, charging or storage, please stop using it and replace it at once. 5. Battery must be removed from the product before discarding it. 6. The product must be powered off before removing the battery. 7. The battery should be safely disposed.
--	--

Component

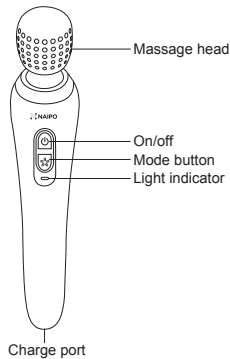

Usage

About charging

1. Before using it for the first time or not using it for a long time, please fully charge it before using it.
2. When the product is in the shutdown state, plug in the USB charging cable to charge it. The indicator light changes from light to dark. When fully charged, the indicator light is always on.
3. When the indicator light flashes during use, it indicates that the battery is insufficient. Please charge it in time.
4. Please unplug the USB charging cable after fully charged.

- ※ Please charge with a mobile phone charger with an output voltage of DC 5V or a computer USB interface for charging.
- ※ This product cannot work when charging. Please disconnect the power supply before use.
- ※ Please keep the battery fully charged every three months if not in use for a long time.

Instruction

On/off :

1. In the shutdown state, press this button to turn on the machine. The indicator light is on and it enters the low-intensity mode by default.
2. In the power-on state, press this button to turn off the machine. The indicator light is off.

Mode button :

1. In the power-on state, press the button and then the indicator light flashes when you switch various working modes. There are 5 working modes (low intensity, medium intensity, high intensity, automatic 1 and automatic 2 for cyclic switching).
2. In the power-on state, long press this button to adjust stepless speed (the speed continues to increase and the users can adjust the speed by themselves). Release it to end the speed adjustment.

When pressing "On/off" button and "Mode" button together :

1. In the power-on state or shutdown state, long press them for 3 seconds to enter the travel lock mode with vibration.
2. In the travel lock mode, long press them for 3 seconds to release the travel lock mode with vibration. You can start using the product then.

Warm tips :

- ※ This massager will be automatically turned on after 15 minutes.
- ※ This product is specially added a travel lock mode to avoid accidental use or battery power consumption, because users may press the buttons by mistake when going out or carrying it. In this mode, all buttons cannot work. After releasing the travel lock mode, the massager can be normally used.

Cleaning and Maintenance

Please turn off the power before cleaning the product.

1. Please clean the product with a slightly damp towel, then dry it with a soft cloth.
2. Please do not use detergent, cleaning brush, gasoline, kerosene, glass cleaner or thinner to clean the product.
3. Please store it in the original box or in a safe, dry and ventilated place after wiping off the dust and dirt.
4. This product is specially designed and manufactured without maintenance or regular lubrication.

Troubleshooting

Please refer to the following tips and possible solutions when you encounter problems during use.

Problems	Possible solutions
Cannot work	<ol style="list-style-type: none"> 1. Ran out of power, please charge it beforeuse. 2. In travel lock mode, please follow the instructions to release travel lock mode.
Stop at work	<ol style="list-style-type: none"> 1. Ran out of power, please charge it before use. 2. The 15-minute timer has expired.
A slight noise at work	It is generated by the mechanical structure, not a malfunction.

If the malfunctions cannot be eliminated according to the above methods, please turn off the power and contact our relevant after-sales service department for help.

About battery recycle

	<ol style="list-style-type: none"> 1. In order to protect the environment, please take out the batteries after the end of their life. Recycle them together, and do not discard them at will. 2. Please do not disassemble the machine at will before the end of its life.
--	--

About battery removal

1. The battery can only be removed when the battery is used up or the massager cannot work at all.
2. Please send the product to the relevant after-sales service department to remove the battery.
3. After removing the battery, please prevent short-circuiting of the positive and negative terminals.

Specification

Product name: Smart handheld massager
 Model: MGPC-806P
 Charging voltage: 5V ---
 Rated voltage: 3.7V ---
 Rated power: 5W

AUKEY INTERNATIONAL LIMITED
 Web: www.naipocare.com
 E-mail: Wholesale@naipocare.com
 Add: No.102, Building P09, China South City, Pinghu Town, Longgang District, Shenzhen, China

