




User Manual

Cordless Handheld Percussion Massager

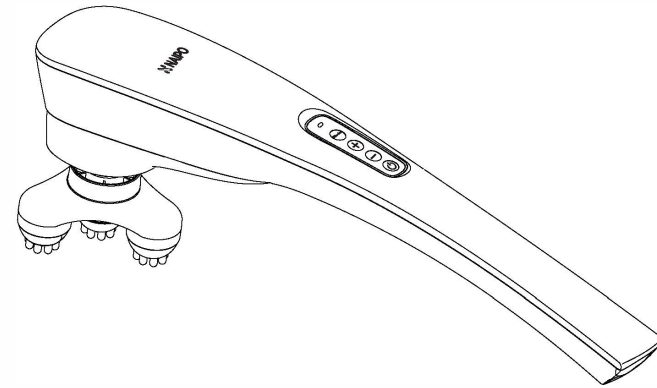

IMPORTANT WARNING:

- The massager is intended for household use only.
- The appliance must only be supplied at SELV.
- This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- Children being supervised not to play with the appliance.
- Cleaning and user maintenance shall not be made by children without supervision.
- Never leave the massager operating or charging unattended.
- This appliance contains batteries that are only replaceable by skilled persons.


- The battery must be removed from the appliance before it is scrapped.
- The appliance must be disconnected from the supply mains when removing the battery.
- The battery is to be disposed of safely.

Removing the built - in rechargeable battery

Remove the built - in battery when disposing of the appliance. Please make sure that the battery is disposed of at an officially designated location if there is one.

This figure must only be used when disposing of the appliance, and must not be used to repair it. If you dismantle the appliance yourself, it may cause malfunction.

- Unplug the adaptor from the appliance before disassembling it.
- Turn on the power and then keep the power on until the battery is completely discharged. Perform steps ① to ② and lift the battery, and then remove it.
- Please take care not to short-circuit the battery. The appliance is only to be used with the power supply unit provided with the appliance.

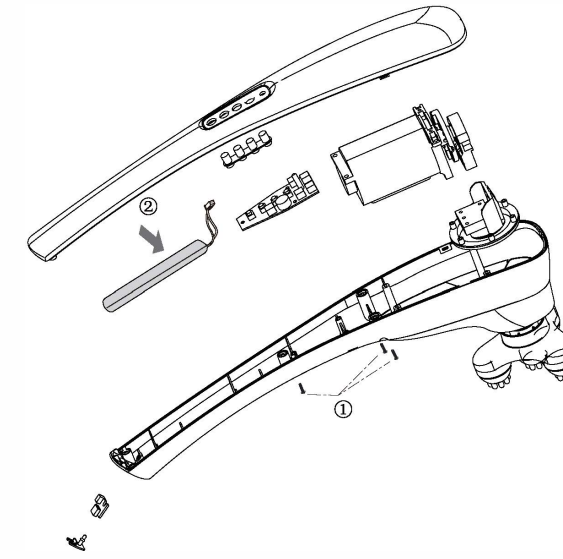

- 01 -

Introduction

Thank you for purchasing the Naipo Cordless Handheld Percussion Massager. Please read this user manual carefully and keep it in a safe place for future reference. If you need any assistance, please contact our support team with your product model number and Amazon order number.

Included

- 1 x Handheld Massager
- 1 x User Manual
- 5 x Interchangeable Attachments
- 1 x power Cable

Product Specs

| | |
|-------------------|---|
| Model | XF-5120 |
| Name | Naipo Cordless Handheld Percussion Massager |
| Power Input | DC 5V |
| Power consumption | 5 W |
| Default time | 20 minutes |
| Li-ion battery | 2200 mAh |
| Charging Time | 2. Hours |
| Usable Time | 1.5 hours |

Important Safeguards

1. Please pull out the plug from the socket when cleaning, maintenance or moving.
2. Please turn off the device if there is something unusual occurs.
3. Please do not use wet hand to pull out or insert the power cord.
4. Please do not use it in humid places such as the bathroom.


- 02 -

5. Please do not destroy, bend or handle the power line, or it may cause fire.
6. Please clear the dust of adaptor regularly.
7. Please avoid expose to sunlight for long time.
8. Please do not modify the product casually.
9. Please do not use it while sleeping.
10. Please do not wash the product by water.
11. The product should be supervised by people if used for long time.
12. Children are prohibited from using this product, please put it in a safe place which is out of their reach.

Questions & Answers


| Problem | Solution |
|---|--|
| The product does not work. | Please check if the power supply is normal. Please reconnect the plug to the mains supply. |
| The product stops working during the use. | The product is overheated because of long-time use. Please disconnect the product from mains power, wait for at least 50 minutes before re-using. Press the power button to restart the product. |
| The product stops working suddenly. | Please check if the power supply is normal. Please reconnect the plug to the mains supply. |

Product Diagram

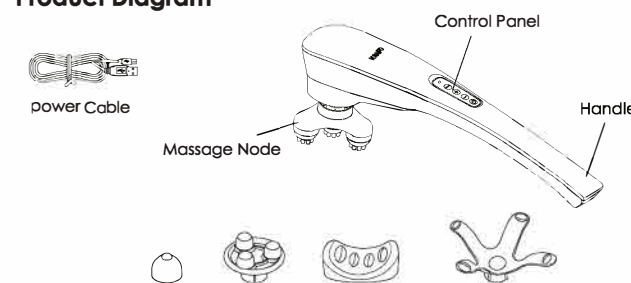

- 03 -

Controller Diagram

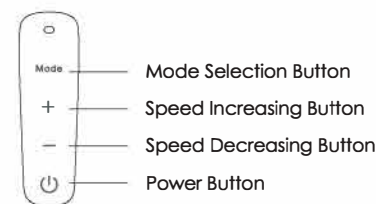

Product Operation Guide

Operation Preparation

1. Select the proper massage node.
2. Connect the charger correctly and plug it into the power supply.

Instructions for Use

Power Button

Press the power button to turn on or turn off the product.

[+]

Press the [+] button to increase the vibration speed [5 optional speed levels].

[-]


Press the [-] button to decrease the vibration speed.

Mode Button

Press the mode selection button to activate automatic massage program [5 optional programs].

- 04 -

Massage Areas

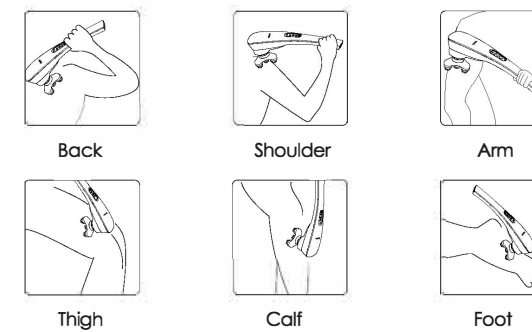

Information on the disposal of Waste Electrical & Electronic Equipment (WEEE) This symbol on the products and accompanying documents means that used electrical and electronic products should not be mixed with general household waste. For proper disposal for treatment, recovery and recycling, please take these products to designated collection points where they will be accepted on a free of charge basis. In some countries you may be able to return your products to your local retailer upon the purchase of a new product. Disposing of this product correctly will help you save valuable resources and prevent any possible effects on human health and the environment, which could otherwise arise from inappropriate waste handling. Please contact your local authority for further details of your nearest collection point for WEEE.

*Please note: Naipo can only provide after sale service for products purchased directly from Naipo.
If you have purchased from a different seller, please contact them directly for service or warranty issues.

- 05 -

Aukey International Limited
Email : Wholesale@naipocare.com
No.102, Building P09, Huanan Town, Pinghu, Longgang District, Shenzhen, 518000 China



Made in China