

材质: 105g铜版纸

尺寸: 105 x 150 mm

印刷: 单黑双面印刷

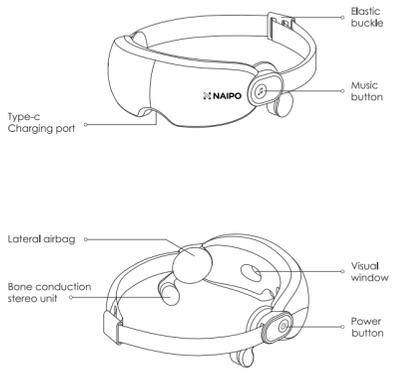
装订方式: 骑马钉

封面封底内页留白, 最

后多出一张为空白页

| | | | |
|---|--|--|--|
| <p>Aukey International Limited Web: www.naipocare.com E-mail: wholesale@naipocare.com Add: 601 6/F, Tung Hip Commercial Building, 248 Des Voeux Road Central, Sheung Wan, Hong Kong</p> <p>CE RoHS    Made in China</p> | <p>NAIPO User Manual</p>  <p>Eye Massager — Model: MGE-2303 —</p> | | <p>Contents</p> <p>Important 01</p> <p>Important Safety Precautions 01</p> <p>Product Diagram 02</p> <p>Usage Method 03</p> <p>Before use 03</p> <p>Indicator light description 03</p> <p>At use 04</p> <p>Maintenance 05</p> <p>How to clean 05</p> <p>Storage 05</p> <p>Error Diagnosis 06</p> <p>Product Specifications 07</p> <p>Warranty & Customer Support 07</p> |
|---|--|--|--|

| |
|---|
| <p>Important</p> <p>Please keep this manual properly for future reference. Before using this appliance, please read all the instructions carefully.</p> <p>Important Safety Precautions</p> <ul style="list-style-type: none">Do not wear accessories around your eyes when using this product.For those with symptoms of eye disease, patients with major eye diseases such as eye surgery, glaucoma, cataract, and retinal detachment, please DO NOT use this product.For those who have medical electronic devices implants, such as cardiac pacemakers, or patients with malignant tumors, heart disease, or severe hypertension, please consult your physician before using this product.Children and mentally/physically disabled can only use the product with special supervision of a responsible adult.If you experience bodily pain and discomfort during use, please stop using it immediately and consult your doctor.This product is not a toy and can not be given to children for fun.This product is for single-person use. Do not operate simultaneously on more than 1 person each time.This product is only for home use. It is not a piece of medical equipment and is not intended for self-diagnosis, treatment, or other commercial purposes.Please charge the product using only accompanied type-c cable and please be careful while you charge.Please cut off the power when the product exhibit malfunctioning. Do not disassemble or fix the product on your own.Do not use the product if you have a fever or body temperature higher than 38 C.Please do not disassemble the product on your own to prevent product malfunctioning.Do not use a power source other than the rated power to prevent fire, physical injuries, or product damages.Do not damage, destroy, over-bend, pull, twist, or tie the type-c cable.Do not expose the product in the environment of high-temperature or hot water.The product should not be discarded with regular home disposal. Please follow the local guideline for disposing of the product.The direction of connecting the product or the type-c cable depends on the product model. Please do not force to connect to the product in case of damaged or shorted connected that may lead to fire or product malfunctions. |
|---|

| |
|--|
| <p>Product Diagram</p>  |
|--|

| <p>Usage Method</p> <p>Before use</p> <ol style="list-style-type: none">Please check the integrity of the machine and accessoriesPlease remove the accessories around the eyes during use to avoid affecting the massage area.Adjust the elastic headband to the size that fits your head shape.Wear the bone conduction stereo unit approximately 0.5-1.5cm in front of your ear. Adjust the different angles of the headphone stand to suit your individual situation for optimal sound.  <p>Indicator light description</p> <table border="1"><thead><tr><th colspan="4">Battery Indicator light</th></tr><tr><th>Color</th><th>Status</th><th colspan="2">Description</th></tr></thead><tbody><tr><td>Orange</td><td>on light</td><td colspan="2">Heating on</td></tr><tr><td>Green</td><td>on light</td><td colspan="2">Heating off</td></tr><tr><td>Red</td><td>breathing light</td><td colspan="2">0%-20%</td></tr><tr><td>Green</td><td>breathing light</td><td colspan="2">20%-99%</td></tr><tr><td>Green</td><td>on light</td><td colspan="2">Fully charged</td></tr><tr><td>Red</td><td>flash light</td><td colspan="2">Low battery</td></tr></tbody></table> | Battery Indicator light | | | | Color | Status | Description | | Orange | on light | Heating on | | Green | on light | Heating off | | Red | breathing light | 0%-20% | | Green | breathing light | 20%-99% | | Green | on light | Fully charged | | Red | flash light | Low battery | |
|---|-------------------------|---------------|--|--|-------|--------|-------------|--|--------|----------|------------|--|-------|----------|-------------|--|-----|-----------------|--------|--|-------|-----------------|---------|--|-------|----------|---------------|--|-----|-------------|-------------|--|
| Battery Indicator light | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | |
| Color | Status | Description | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | |
| Orange | on light | Heating on | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | |
| Green | on light | Heating off | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | |
| Red | breathing light | 0%-20% | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | |
| Green | breathing light | 20%-99% | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | |
| Green | on light | Fully charged | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | |
| Red | flash light | Low battery | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | |

| Music indicator light | | | |
|----------------------------------|-------|-------------|----------------------------------|
| | Color | Status | Description |
| Normal working status | Blue | on light | Bluetooth connected |
| | Cyan | on light | Bluetooth not connected |
| Search for new Bluetooth devices | Blue | flash light | Search for new Bluetooth devices |

At use

| Buttons | |
|---|--|
|  | 1. Long press the power button 2 seconds, ON/OFF |
| | 2. Double-click the power button (press twice within 1 seconds) to turn on/off the heating. |
| | 3. Default power on into the "Energy Mode", short press the power button to switch the massage mode, in order of "Energy Mode", "Relax Mode", "Relieve Mode", "Sleep Mode". |
|  | 1. The first time to connect Bluetooth, long press the music button 2 seconds, search for new Bluetooth devices, the Bluetooth name for this product is "NAIPO". If need to disconnect current Bluetooth device, please long press the music button for 2 seconds. |
| | 2. Double-click the music button (press twice within 1 seconds) to turn on/off the music. |
| | 3. Short press the music button to switch to the next music. |

*Without connecting Bluetooth status, massage mode runs for 15 minutes and ends, connecting Bluetooth status, massage mode runs for 15 minutes and ends, music function runs for 30 minutes and ends

| |
|--|
| <p>Maintenance</p> <p>How to clean</p> <ol style="list-style-type: none">The product shell can be cleaned with a clean cloth;Stains can be removed by rubbing alcohol on a clean rag;For deep dirt, wet the cloth with neutral detergent and gently wipe it off.Do not pour water to the productLeave it to dry naturally. <p>Storage</p> <p>Please place the product inside of the storage bag with its original package and store it in a safe place.</p> <p>Do not put heavy weights on top of the product.</p> <p>Please take moisture, rodent and pollution prevention measures before storing the product.</p> <p>For product in more than 3 months of storage, please recharge before use.</p> <p>Do not store the product in places with direct sunlight, or corrosive gas.</p> |
|--|

| <p>Error Diagnosis</p> <table border="1"><thead><tr><th>Error</th><th>Reason</th><th>Solution</th></tr></thead><tbody><tr><td>Can not turn on.</td><td>Insufficient battery.</td><td>Please recharge.</td></tr><tr><td>Air-pressure strength is weak.</td><td>The strap is too loose.</td><td>Adjust the length of the strap, make the massager comfortable to wear</td></tr><tr><td rowspan="2">Stopped working during operation.</td><td>Massage timer is off.</td><td>Press the Power button to reactivate the product.</td></tr><tr><td>Battery is ran out.</td><td>Recharge the battery.</td></tr><tr><td>Cannot connect Bluetooth</td><td>Mistakenly touched the music button. Bluetooth is turned off</td><td>Long press the power button for 2 seconds to turn off, and the Bluetooth function will be restored after restarting.</td></tr><tr><td>Bluetooth playback music sound is low</td><td>Mobile phone playback volume is low</td><td>Adjusting the phone volume on your phone</td></tr></tbody></table> | Error | Reason | Solution | Can not turn on. | Insufficient battery. | Please recharge. | Air-pressure strength is weak. | The strap is too loose. | Adjust the length of the strap, make the massager comfortable to wear | Stopped working during operation. | Massage timer is off. | Press the Power button to reactivate the product. | Battery is ran out. | Recharge the battery. | Cannot connect Bluetooth | Mistakenly touched the music button. Bluetooth is turned off | Long press the power button for 2 seconds to turn off, and the Bluetooth function will be restored after restarting. | Bluetooth playback music sound is low | Mobile phone playback volume is low | Adjusting the phone volume on your phone |
|--|--|--|----------|------------------|-----------------------|------------------|--------------------------------|-------------------------|---|-----------------------------------|-----------------------|---|---------------------|-----------------------|--------------------------|--|--|---------------------------------------|-------------------------------------|--|
| Error | Reason | Solution | | | | | | | | | | | | | | | | | | |
| Can not turn on. | Insufficient battery. | Please recharge. | | | | | | | | | | | | | | | | | | |
| Air-pressure strength is weak. | The strap is too loose. | Adjust the length of the strap, make the massager comfortable to wear | | | | | | | | | | | | | | | | | | |
| Stopped working during operation. | Massage timer is off. | Press the Power button to reactivate the product. | | | | | | | | | | | | | | | | | | |
| | Battery is ran out. | Recharge the battery. | | | | | | | | | | | | | | | | | | |
| Cannot connect Bluetooth | Mistakenly touched the music button. Bluetooth is turned off | Long press the power button for 2 seconds to turn off, and the Bluetooth function will be restored after restarting. | | | | | | | | | | | | | | | | | | |
| Bluetooth playback music sound is low | Mobile phone playback volume is low | Adjusting the phone volume on your phone | | | | | | | | | | | | | | | | | | |

| | | | | | | | | | | | | | | |
|---|--------------|--------------|---------------|----------|---------------|----|------------------|----|------------------|---------|--------------|--------------|------------|------|
| <p>Product Specifications</p> <table border="1"><tbody><tr><td>Product name</td><td>Eye Massager</td></tr><tr><td>Product model</td><td>MGE-2303</td></tr><tr><td>Product power</td><td>4W</td></tr><tr><td>charging voltage</td><td>5V</td></tr><tr><td>battery capacity</td><td>1200mAh</td></tr><tr><td>Product size</td><td>217×114×80mm</td></tr><tr><td>Net weight</td><td>360g</td></tr></tbody></table> <p>Warranty & Customer Support</p> <p>If you need help regarding product issues, after-sales service, or other relevant problems, please contact us via the email on the cover.</p> <p>*Note: Products purchased directly from NAIPO benefit from our after-sales service. If the purchases are made from local retail stores, we recommend you contact them for warranty issues.</p> | Product name | Eye Massager | Product model | MGE-2303 | Product power | 4W | charging voltage | 5V | battery capacity | 1200mAh | Product size | 217×114×80mm | Net weight | 360g |
| Product name | Eye Massager | | | | | | | | | | | | | |
| Product model | MGE-2303 | | | | | | | | | | | | | |
| Product power | 4W | | | | | | | | | | | | | |
| charging voltage | 5V | | | | | | | | | | | | | |
| battery capacity | 1200mAh | | | | | | | | | | | | | |
| Product size | 217×114×80mm | | | | | | | | | | | | | |
| Net weight | 360g | | | | | | | | | | | | | |

| |
|--|
| |
|--|